




ken-a-vision
KNOWLEDGE THROUGH VISION

Video Flex[®]

7600

The Video Flex 7600 introduces a direct computer interface. Simply connect to the USB port of any PC to transmit or capture images. The 7600 can simultaneously display to a computer and a TV monitor. This means that a video signal can be received and "frozen" in a computer while live action continues on the monitor.

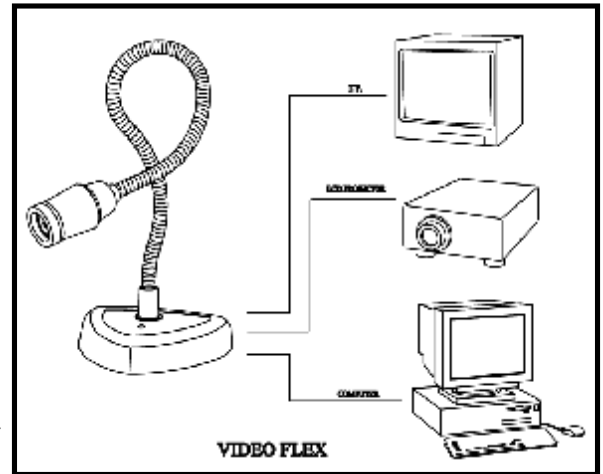

Additional Features and Benefits

- **USB, Composite, & S-VHS** For greater range of applications; simultaneous hook up
- **Triangular Base** Slip and Tip proof
- **500 Lines of Resolution** For crisp clear pictures
- **Poly-Carbonate Build** Allows for rugged conditions and environment
- **C-mount Lens** Interchangeable lens system accommodates all c-mount lenses.
- **Quick Focus™** 1/4 turn moves from 6mm to infinity
- **Ball and Socket Head** For greater range of motion and easier positioning
- **Microphone** For recording or amplifying
- **Discovery Scope Kit** To bring nature into the classroom
- **Adjustable Contrast** Varies the contrast of the picture
- **Positive / Negative Switch** Similar to Darkfield Control on microscopes; this reverses the lighted images
- **Fluorescent Control** For Fluorescent conditions
- **KensingtonSlot** Provides anti-theft control

7600 SPECIFICATIONS

Weight	5.95 lbs / 2.70 Kg
Base Material	Polycarbonate
Base Dimension	7"x7"x7" / 18x18x18 cm
Height	30" / 76.2 cm
Head Type	Ball and Socket
Lens	8 mm
Focus	¼ turn "Quick Focus"
Focal distance	¼" / 6mm to infinity
Light Sensitivity	1.5 Lux
Resolution	500 Lines
Magnification	50:1 and greater

Exposure	Adjustable Contrast
Output Signal	USB (cable and software included)
	Composite and S-VHS
Eyepiece Adapters	28mm & Built-in 34.5 mm
Advanced Features	Contrast Control Buttons
	Positive/Negative Switch
	Fluorescent Control Switch
Options Included	Discovery Scope Kit & Case
Cable	12' A/V
Warranty	5 Years

Specifications subject to change without notice



ken-a-vision
KNOWLEDGE THROUGH VISION

(816) 353-4787 Fax (816) 358-5072

email: info@ken-a-vision.com